
	<b>TITLE: DEVICE MASTER RECORD</b>	
	<b>QUALITY FILE</b>	<b>Ref. no.: QF 04-5</b>
	<b>Revision: 1</b>	<b>Page 1 of 1</b>
	<b>Issue date: 04.04.2022</b>	<b>Effective date: 04.04.2022</b>

## SUCTION SET

REF. NO. 1148 3532 80 ML

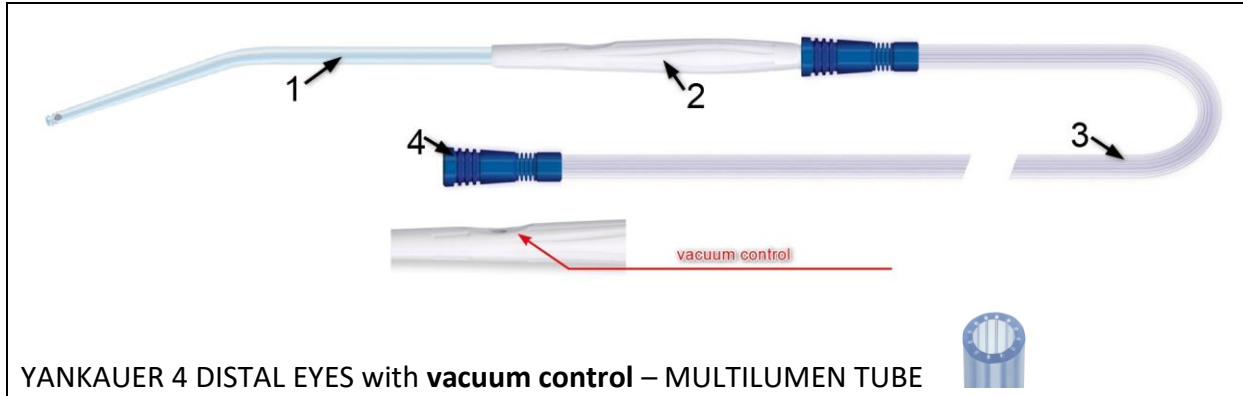

### COMPONENT DESCRIPTION

### COMPONENT NO.

#### 1. SUCTION PIECE

rigid tube ID/OD 3.0/4.0, 186 mm, 4 distal eyes, bent  
material: PVC

0.48022040

#### 2. HANDLE SHORT

handle white ID 6.5 mm, 160 mm length, with **vacuum control**  
material: Styrolux (SBC)

4.20016065

#### 3. TUBE

Multi-Lumen transparent tube ID/OD 5.6/8.0 mm,  
material: PVC 73°ShoreA, DEHP – FREE,  
• 1148 3532 80 ML – 3400 mm length

0.56340080

#### 3. PROXIMAL END

funnel 8.0 PVC  
material: PVC 50°ShoreA, DEHP – FREE, blue

0.00000080

### PACKAGING MATERIAL

Innerbag 180 x 380 mm		2.18038004
Peelpack (packing machine) 200 x 350 mm		5.20035410
Outerbag 310 x 600 mm	<b>25 pcs/bag</b>	2.31060010
Box 386 x 386 x 260 mm	<b>50 pcs/box</b>	6.38638626
Label on box 85 x 131 mm		7.00085131

Work instruction	QI 01
Packaging instruction	QI 04

### Comments:

	Glued	Not glued
Suction piece and handle	X	
<b>Suction piece and handle short put into innerbag.</b>		