



CoperSurgical


Lone Star Reusable Retractor Rings

For simpler, better managed surgeries.

For a wide range of surgical procedures

Benefits

- Excellent visibility
- Adapts and repositions quickly
- Adjustable tension for optimal exposure
- Frees up OR staff

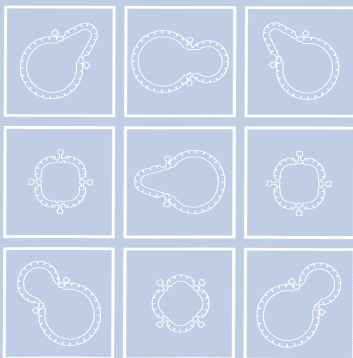

Maximum exposure.

Minimum effort.

Resilient and reusable.


Lone Star Reusable Retractor Rings

Flexible. Durable. Remarkable!

So simple - and yet so capable. These Reusable Retractor Rings have the qualities you need for optimum surgical results. They dramatically improve visibility in the surgical field, by eliminating the need for manual retractors and the personnel to hold them. The articulated design adapts to any region of the anatomy. Adjustable hinges and a system of Elastic Stay hooks allow easy positioning and repositioning throughout the procedure, without placing undue pressure on delicate tissues.

Resilient aluminum and durable stainless steel construction stand up to repeated use. Every detail supports more efficient surgeries from start to finish.

Ease and Precision in 3 Simple Steps

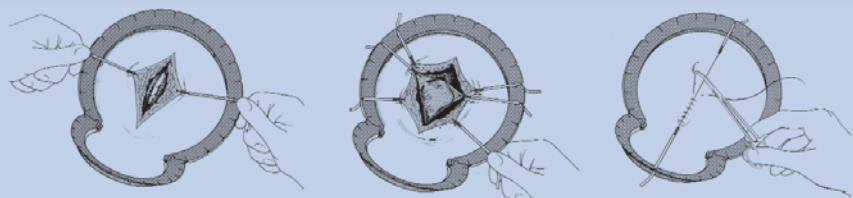

1. Place Ring around initial incision. Insert Elastic Stay hooks into tissue and retract in directly opposing quadrants. Slots on ring allow easy placement for proper tension.
2. Reposition or add hooks throughout procedure to maintain balanced retraction.
3. When closing, place stays along axis of the incision to create tension to facilitate precise staple or suture placement.

Reusable Retractor Rings Specifications

Stainless Steel Retractor Rings

Packaged singly, ready for sterilization.

2200 27.5 cm x 18.4 cm (795 grams)

2201 16.4 cm x 11.6 cm (340 grams)

Aluminum Retractor Rings

Packaged singly, ready for sterilization.

1103 11.6 cm (96 grams)

1104 32.2 cm x 18.4 cm (236 grams)


1107 14.5 cm x 14.5 cm (206 grams)

1110 24.3 cm x 24.3 cm (380 grams)



1107

1104

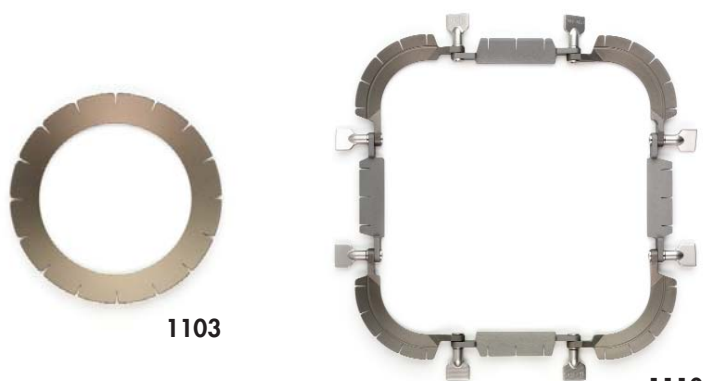

1103

1110



2201

2200

CooperSurgical

www.coopersurgical.com

To find out more about the Lone Star Reusable Retractor Rings, or to place an order call 800.243.2974.

Lone Star Retractor System™ is a trademark of CooperSurgical, Inc.

95 Corporate Drive
Trumbull, CT 06611
800.243.2974
203.601.5200